
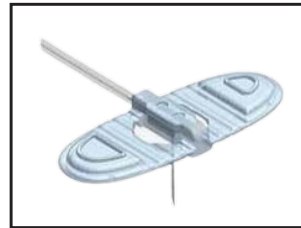
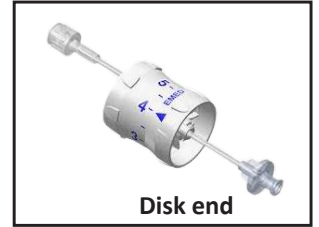


SCIG Preparation Guide





Preparing the Tubing

1. Remove the cap from the disc end of the flow control tubing and connect the filled medication syringe.
2. Attach needle set to other end of flow control tubing.
3. Prime the infusion set before connecting the needle sets if instructed to do so.
4. Do not prime subcutaneous needles as a SCIG requires a dry insertion technique to avoid skin irritation.
5. Prepare your infusion site(s). See SCIG Site Preparation Guide.

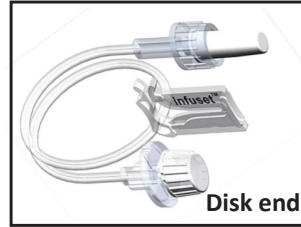

Preparing the Pump

1. Unwind blue end.
2. Insert syringe and twist.
3. Flange locks into slot.
4. Twist blue knob until flush with white portion.
5. Unclamp tubing.

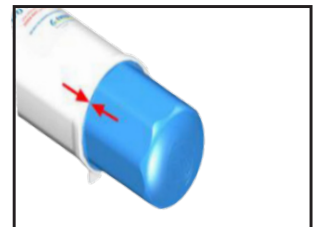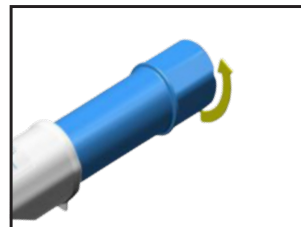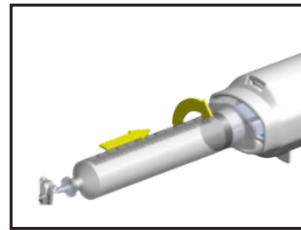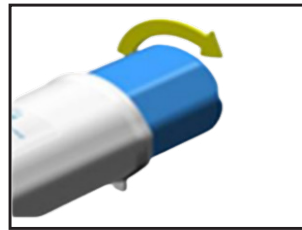

If using more than one syringe, repeat pump directions 1-5 and switch syringes.

Infusion Complete

1. Remove syringe from the pump by unwinding the blue knob.
2. While holding the needle set in place, gently remove the dressing, fold up the wings, and pull the needle(s) straight out and place them in sharps container.
3. Press a small piece of gauze on the needle site and apply tape or adhesive bandage.
4. Record the infusion in your infusion log.